

Side Guard Installation Guide

Have a question? Call 617-356-8106 or email support@gogrit.us.

The Side Guards will come with both magnets and screws. Most folks prefer to use the magnets only (for easier removal), but inserting just one screw will secure the Side Guard to the seatpan. Note: Both Side Guards are identical: There is no “left” or “right.”

To install the Side Guards:

1. Position the Side Guard as shown in Figure 1. You may find it easier to remove the seat cushion to ensure proper alignment.
2. The magnets will secure the Side Guard to the sidewall of the seatpan. If you'd like to install a screw as well, line up the Side Guard with the corresponding hole on the seatpan (see Figure 2).
3. Slide the washer over the threaded portion of the screw. Position the screw so that the bolt head is on the outside of the seatpan. Tighten the screw with a Phillips head screwdriver.
4. Repeat for the other Side Guard.



Figure 1: Shown without seat cushion.

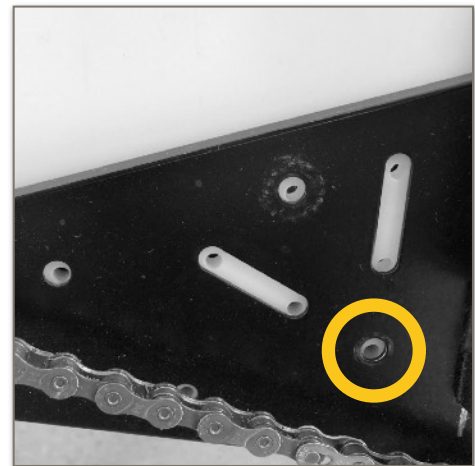


Figure 2: Insert the screw through this hole.



Figure 3: Finished look at an installed Side Guard.